



# Histone H3 (Mono Methyl Lys37) rabbit pAb

Cat No.:ES1080

For research use only

## Overview

<b>Product Name</b>	Histone H3 (Mono Methyl Lys37) rabbit pAb
<b>Host species</b>	Rabbit
<b>Applications</b>	WB;IHC;IF;ELISA
<b>Species</b>	Human;Mouse;Rat
<b>Cross-Reactivity</b>	
<b>Recommended dilutions</b>	Western Blot: 1/500 - 1/2000. IHC-p: 1:100-300 ELISA: 1/20000. Not yet tested in other applications.
<b>Immunogen</b>	Synthesized peptide derived from human Histone H3 around the mono-methylation site of K37.
<b>Specificity</b>	Mono-Methyl-Histone H3 (K37) Polyclonal Antibody detects endogenous levels of Histone H3 protein only when mono-methylated at K37.
<b>Formulation</b>	Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide.
<b>Storage</b>	Store at -20°C. Avoid repeated freeze-thaw cycles.
<b>Protein Name</b>	Histone H3.1/Histone H3.2/Histone H3.3/Histone H3.3C
<b>Gene Name</b>	HIST1H3A/HIST1H3/HIST1H3C/HIST1H3D/HIST1H3E/HIST1H3F/HIST1H3G/HIST1H3H/HIST1H3I/HIST1H3J/HIST2H3A/HIST2H3C/HIST2H3D/H3F3A/H3F3B/H3F3C
<b>Cellular localization</b>	Nucleus. Chromosome.
<b>Purification</b>	The antibody was affinity-purified from rabbit antiserum by





**ation** affinity-chromatography using epitope-specific immunogen.

**Clonality** Polyclonal

**Concentration** 1 mg/ml

**Observed band** 17kD

**Human Gene ID** 8350/8351/8352/8353/8354/8355/8356/8357/8358/8968/126961/33393

**Human Swiss-Prot Number** P68431/Q71DI3/P84243/Q6NXT2


**Alternative Names** H3K37ME1; HIST1H3A; H3FA; HIST1H3B; H3FL; HIST1H3C; H3FC; HIST1H3D; H3FB; HIST1H3E; H3FD; HIST1H3F; H3FI; HIST1H3G; H3FH; HIST1H3H; H3FK; HIST1H3I; H3FF; HIST1H3J; H3FJ; Histone H3.1; Histone H3.1; Histone H3.2; Histone H3/m; Histone H3/o; H3F3A; H3.3A; H3

**Background** Histones are basic nuclear proteins that are responsible for the nucleosome structure of the chromosomal fiber in eukaryotes. This structure consists of approximately 146 bp of DNA wrapped around a nucleosome, an octamer composed of pairs of each of the four core histones (H2A, H2B, H3, and H4). The chromatin fiber is further compacted through the interaction of a linker histone, H1, with the DNA between the nucleosomes to form higher order chromatin structures. This gene is intronless and encodes a replication-dependent histone that is a member of the histone H3 family. Transcripts from this gene lack polyA tails; instead, they contain a palindromic termination element. This gene is found in the large histone gene cluster on chromosome 6p22-p21.3. [provided by RefSeq, Aug 2015],






**ELK Biotechnology**

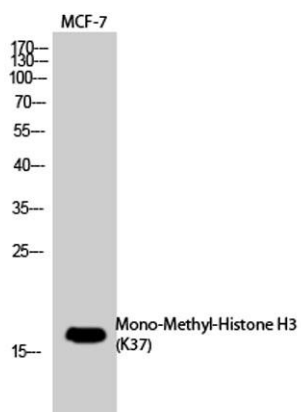


Western Blot analysis of MCF-7, PC-12 cells using Mono-Methyl-Histone H3 (K37) Polyclonal Antibody. Secondary antibody(catalog#:RS0002) was diluted at 1:20000

Immunohistochemical analysis of paraffin-embedded human-breast, antibody was diluted at 1:100

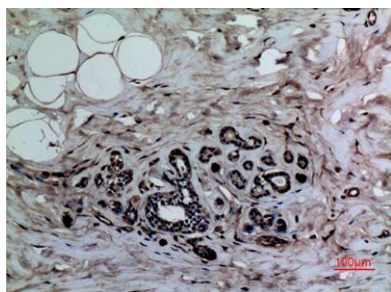

Immunohistochemical analysis of paraffin-embedded human-breast, antibody was diluted at 1:100

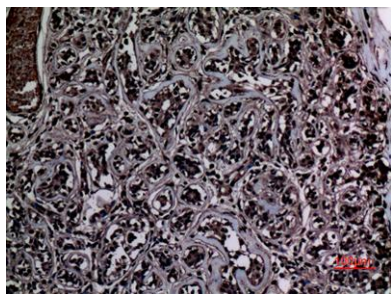

+86-27-59760950

[ELKbio@ELKbiotech.com](mailto:ELKbio@ELKbiotech.com)

[www.elkbiotech.com](http://www.elkbiotech.com)

23-2, No.388 Gaoxin 2nd Road,Wuhan East Lake Hi-tech Development Zone, Hubei , P.R.C