

CXCR4 rabbit pAb

Cat No.: ES20775

For research use only

Overview

Product Name CXCR4 rabbit pAb

Host species Rabbit
Applications WB;IHC;IF

Species Cross-Reactivity Human; Rat; Mouse

Recommended dilutions WB 1:1000-2000, IHC 1:100-200

ImmunogenSynthetic Peptide of CXCR4 AA range: 203-253SpecificityCXCR4 protein(A220) detects endogenous levels of

CXCR4

Formulation Liquid in PBS containing 50% glycerol, 0.5% BSA and

0.02% sodium azide.

Storage Store at -20°C. Avoid repeated freeze-thaw cycles.

Protein Name C-X-C chemokine receptor type 4 (CXC-R4) (CXCR-4)

(FB22) (Fusin) (HM89) (LCR1) (Leukocyte-derived seven transmembrane domain receptor) (LESTR) (NPYRL) (Stromal cell-derived factor 1 receptor)

(SDF-1 r

Gene Name CXCR4

Cell membrane ; Multi-pass membrane protein . Cell

junction. Early endosome. Late endosome.

Lysosome. In unstimulated cells, diffuse pattern on

plasma membrane. On agonist stimulation, colocalizes with ITCH at the plasma membrane where it becomes ubiquitinated. In the presence of antigen, distributes to the immunological synapse forming at the T-cell-APC contact area, where it localizes at the peripheral and distal supramolecular

activation cluster (SMAC).

Purification The antibody was affinity-purified from rabbit

antiserum by affinity-chromatography using

epitope-specific immunogen.

ClonalityPolyclonalConcentration1 mg/mlObserved band50kD



+86-27-59760950 ELKbio@ELKbiotech.com www.elkbiotech.com



Human Gene ID Human Swiss-Prot Number Alternative Names 7852 P61073


CXCR4; C-X-C chemokine receptor type 4; CXC-R4;

CXCR-4; FB22; Fusin; HM89; LCR1;


Leukocyte-derived seven transmembrane domain receptor; LESTR; NPYRL; Stromal cell-derived factor 1 receptor; SDF-1 receptor; CD antigen CD184 C-X-C motif chemokine receptor 4(CXCR4) Homo This gene encodes a CXC chemokine sapiens receptor specific for stromal cell-derived factor-1. The protein has 7 transmembrane regions and is located on the cell surface. It acts with the CD4 protein to support HIV entry into cells and is also highly expressed in breast cancer cells. Mutations in this gene have been associated with WHIM (warts, hypogammaglobulinemia, infections, and myelokathexis) syndrome. Alternate transcriptional splice variants, encoding different isoforms, have been characterized. [provided by RefSeq, Jul 2008],

Background

Fat Hela brain Western blot analysis of 1) Rat Brain Tissue, 2)Hela, with CXCR4 Rabbit pAb diluted at 1:2,000.

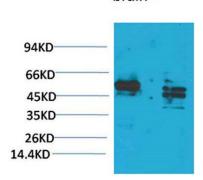


Immunohistochemical analysis of paraffin-embedded Human Breast Carcinoma using CXCR4Rabbit pAb diluted at 1:200.

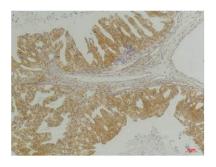





Immunohistochemical analysis of paraffin-embedded Human Tonsil Tissue using CXCR4Rabbit pAb diluted at 1:200.

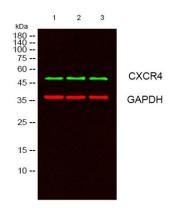

Western blot analysis of lysates from 1) Rat Brain Tissue, 2)Hela, 3) HEK293 cells, (Green) primary antibody was diluted at 1:1000, 4°over night, secondary antibody(cat:RS23920)was diluted at 1:10000, 37° 1hour. (Red) GAPDH Monoclonal Antibody(2B8) (cat:YM3029) antibody was diluted at 1:5000 as loading control, 4° over night, secondary antibody(cat:RS23710)was diluted at 1:10000, 37° 1hour.

