

## Ksr-1 rabbit pAb

Cat No.: ES7889

For research use only

## Overview

Product Name Ksr-1 rabbit pAb

Host species Rabbit

Applications WB;IHC;IF;ELISA Species Cross-Reactivity Human;Mouse

**Recommended dilutions** Western Blot: 1/500 - 1/2000.

Immunohistochemistry: 1/100 - 1/300.

Immunofluorescence: 1/200 - 1/1000. ELISA: 1/20000. Not yet tested in other applications.

Immunogen The antiserum was produced against synthesized

peptide derived from human KSR. AA range:358-407

**Specificity** Ksr-1 Polyclonal Antibody detects endogenous levels

of Ksr-1 protein.

Formulation Liquid in PBS containing 50% glycerol, 0.5% BSA and

0.02% sodium azide.

**Storage** Store at -20°C. Avoid repeated freeze-thaw cycles.

Protein Name Kinase suppressor of Ras 1

Gene Name KSR1

Cellular localization Cytoplasm . Membrane ; Peripheral membrane

protein . Cell membrane ; Peripheral membrane

protein. Cell projection, ruffle membrane.

Endoplasmic reticulum membrane . In unstimulated cells, where the phosphorylated form is bound to a 14-3-3 protein, sequestration in the cytoplasm occurs. Following growth factor treatment, the protein is free for membrane translocation, and it moves from the cytoplasm to the cell periphery. .

Purification The antibody was affinity-purified from rabbit

antiserum by affinity-chromatography using

epitope-specific immunogen.

ClonalityPolyclonalConcentration1 mg/mlObserved band115kDHuman Gene ID8844



+86-27-59760950 ELKbio@ELKbiotech.com

www.elkbiotech.com




Human Swiss-Prot Number Alternative Names Background Q8IVT5


KSR1; KSR; Kinase suppressor of Ras 1 caution: The sequence shown here is derived from an Ensembl automatic analysis pipeline and should be considered as preliminary data., function: Location-regulated scaffolding protein connecting MEK to RAF. Promotes MEK and RAF phosphorylation and activity through assembly of an activated signaling complex. By itself, it has no demonstrated kinase activity., PTM: Phosphorylated on Ser-309 and, to a higher extent, on Ser-404 by MARK3. Dephosphorylated on Ser-404 by PPP2CA. In resting cells, phosphorylated KSR1 is cytoplasmic and in stimulated cells, dephosphorylated KSR1 is membrane-associated., similarity: Belongs to the protein kinase superfamily. TKL Ser/Thr protein kinase family., similarity: Contains 1 phorbol-ester/DAG-type zinc finger., similarity: Contains 1 protein kinase domain., subcellular location: In unstimulated cells, where the phosphorylated form is bound to a 14-3-3 protein, sequestration in the cytoplasm occurs. Following growth factor treatment, the protein is free for membrane translocation, and it moves from the cytoplasm to the cell periphery., subunit: Interacts with HSPCA/HSP90, YWHAB/14-3-3, CDC37, MAP2K/MEK, MARK3, PPP2R1A and PPP2CA. Also interacts with RAF and MAPK/ERK, in a Ras-dependent manner (By similarity). The binding of 14-3-3 proteins to phosphorylated KSR prevents the membrane localization.,






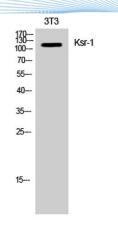


Western Blot analysis of NIH-3T3 cells using Ksr-1 Polyclonal Antibody

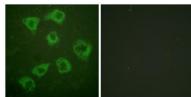

Immunofluorescence analysis of HUVEC cells, using KSR Antibody. The picture on the right is blocked with the synthesized peptide.

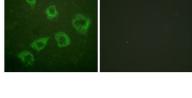


Immunohistochemistry analysis of paraffin-embedded human brain tissue, using KSR Antibody. The picture on the right is blocked with the synthesized peptide.

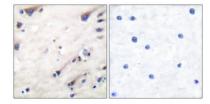





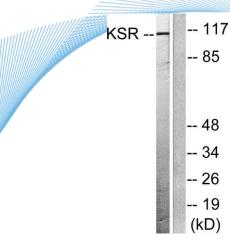

Western blot analysis of lysates from NIH/3T3 cells, treated with PDGF 50ng/ml 20', using KSR Antibody. The lane on the right is blocked with the synthesized peptide.

