



Recombinant Human CFHR3 (C-6His)

Catalog #	EPT083
Expression Host	Human Cells
DESCRIPTION	Recombinant Human Complement Factor H-related Protein 3 is produced by our Mammalian expression system and the target gene encoding Gln19-Glu330 is expressed with a 6His tag at the C-terminus.
Accession	Q02985
Synonyms	Complement factor H-related protein 3; FHR-3; DOWN16; H factor-like protein 3; CFHR3; CFHL3
Mol Mass	36.4 KDa
AP Mol Mass	43-56 KDa, reducing conditions
Purity	Greater than 95% as determined by reducing SDS-PAGE.
Endotoxin	Less than 0.1 ng/μg (1 EU/μg) as determined by LAL test.
FORMULATION	Lyophilized from a 0.2 μm filtered solution of PBS, pH 7.4.
RECONSTITUTION	Always centrifuge tubes before opening. Do not mix by





vortex or pipetting.

It is not recommended to reconstitute to a concentration less than 100µg/ml.

Dissolve the lyophilized protein in distilled water.

Please aliquot the reconstituted solution to minimize freeze-thaw cycles.

SHIPPING

The product is shipped at ambient temperature.

Upon receipt, store it immediately at the temperature listed below.

STORAGE

Lyophilized protein should be stored at $< -20^{\circ}\text{C}$, though stable at room temperature for 3 weeks.

Reconstituted protein solution can be stored at $4-7^{\circ}\text{C}$ for 2-7 days.

Aliquots of reconstituted samples are stable at $< -20^{\circ}\text{C}$ for 3 months.

BACKGROUND

Complement factor H-related protein 3 is a protein that in humans is encoded by the CFHR3 gene. It is a secreted protein that belongs to the complement factor H protein family. Members of the H-related protein family are exclusively composed of individually folded protein domains. CFHR3 is synthesized as a 330 amino acid precursor that contains an 18 amino acid






ELK Biotechnology

signal peptide and a 312 amino acid mature chain.

Human CF3R3 protein might be involved in complement regulation.

SDS-PAGE

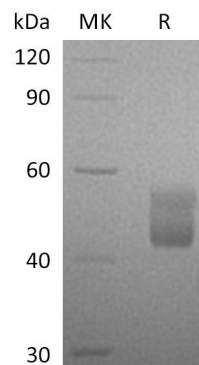

+86-27-59760950

ELKbio@ELKbiotech.com

www.elkbiotech.com

23-2, No.388 Gaoxin 2nd Road, Wuhan East Lake Hi-tech Development Zone, Hubei, P.R.C