



Recombinant Human Mindin (C-6His)

Catalog #	EPT104
Expression Host	Human Cells
DESCRIPTION	Recombinant Human Spondin2 is produced by our Mammalian expression system and the target gene encoding Gln27-Val331 is expressed with a 6His tag at the C-terminus.
Accession	AAH02707.1
Synonyms	Spondin-2; Differentially expressed in cancerous and non-cancerous lung cells 1; DIL-1; Mindin; SPON2
Mol Mass	34.4 KDa
AP Mol Mass	38-42 KDa, reducing conditions
Purity	Greater than 95% as determined by reducing SDS-PAGE.
Endotoxin	Less than 0.1 ng/μg (1 EU/μg) as determined by LAL test.
FORMULATION	Lyophilized from a 0.2 μm filtered solution of 20mM PB, 150mM NaCl, 1mM EDTA, pH 7.4.
RECONSTITUTION	Always centrifuge tubes before opening. Do not mix by





vortex or pipetting.

It is not recommended to reconstitute to a concentration less than 100µg/ml.

Dissolve the lyophilized protein in distilled water.

Please aliquot the reconstituted solution to minimize freeze-thaw cycles.

SHIPPING

The product is shipped at ambient temperature.

Upon receipt, store it immediately at the temperature listed below.

STORAGE

Lyophilized protein should be stored at $< -20^{\circ}\text{C}$, though stable at room temperature for 3 weeks.

Reconstituted protein solution can be stored at $4-7^{\circ}\text{C}$ for 2-7 days.

Aliquots of reconstituted samples are stable at $< -20^{\circ}\text{C}$ for 3 months.


BACKGROUND

Spondin-2, also referred to as mindin, belongs to the F-spondin family of secreted extracellular matrix proteins. Spondins are characterised by the presence of F-spondin domains 1 and 2 (FS1 and FS2) at the N-terminus and a thrombospondin-type 1 repeat (TSR1) domain at the C-terminus. Spondin-2 functions as a pattern-recognition molecule for bacterial and





viral pathogens and as an integrin ligand for inflammatory cell recruitment and T cell priming. In addition to its roles in promoting neuron outgrowth and inhibiting both cancer and angiogenesis, Spondin-2 plays an important role in the initiation of the immune response and is involved in inflammatory processes.

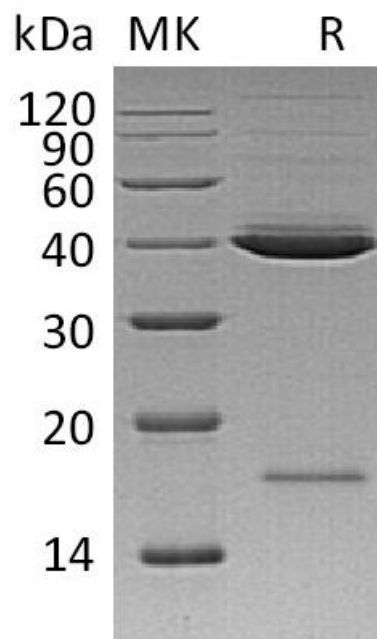

SDS-PAGE

