

Recombinant Mouse SLAMF2 (C-6His)

Catalog # EPT206

Expression Host Human Cells

DESCRIPTION Recombinant Mouse SLAM Family Member 2 is

produced by our Mammalian expression system and

the target gene encoding Phe23-Arg216 is expressed

with a 6His tag at the C-terminus.

Accession Q18PH8

Synonyms CD48; BLAST1; BCM1; SLAMF2;SLAM Family Member 2

Mol Mass 23.1 KDa

AP Mol Mass 45-50 KDa, reducing conditions

Purity Greater than 95% as determined by reducing

SDS-PAGE.

Endotoxin Less than 0.1 ng/μg (1 EU/μg) as determined by LAL

test.

FORMULATION Lyophilized from a 0.2 µm filtered solution of PBS, pH

7.4.

RECONSTITUTION Always centrifuge tubes before opening.Do not mix by

vortex or pipetting.



www.elkbiotech.com



It is not recommended to reconstitute to a concentration less than 100µg/ml.

Dissolve the lyophilized protein in distilled water.

Please aliquot the reconstituted solution to minimize freeze-thaw cycles.

SHIPPING

The product is shipped at ambient temperature.

Upon receipt, store it immediately at the temperature listed below.

STORAGE

Lyophilized protein should be stored at < -20 ° C, though stable at room temperature for 3 weeks.

Reconstituted protein solution can be stored at 4-7°C for 2-7 days.

Aliquots of reconstituted samples are stable at < -20° C for 3 months.


BACKGROUND

CD48 is a GPI-linked protein in the CD2 family of immunoglobulin superfamily proteins. CD48 is expressed on most lineage-committed hematopoietic cells but not on hematopoietic stem cells or multipotent hematopoietic progenitors. Among dendritic cells (DC), CD48 is selectively expressed on circulating myeloid DC and resident bone marrow and thymus DC. CD2, 2B4, and heparan sulfate function as





CD48 ligands. CD48 is competent to transduce signals and can also trigger signaling through CD2 or 2B4. CD48 expressed on NK cells is coactivating, whereas CD48 expressed on other cell types inhibits NK cell activation.

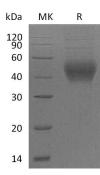

SDS-PAGE

