



Recombinant Human PDYN (C-6His)

Catalog #	EPT300
Expression Host	Human Cells
DESCRIPTION	Recombinant Human Proenkephalin-B is produced by our Mammalian expression system and the target gene encoding Asp21-Ala254 is expressed with a 6His tag at the C-terminus.
Accession	P01213
Synonyms	ADCA; PENKB; SCA23; PDYN; Dynorphin B; Big Dyn
Mol Mass	27.3 KDa
AP Mol Mass	35 KDa, reducing conditions
Purity	Greater than 95% as determined by reducing SDS-PAGE.
Endotoxin	Less than 0.1 ng/μg (1 EU/μg) as determined by LAL test.
FORMULATION	Lyophilized from a 0.2 μm filtered solution of 4mM HCl.
RECONSTITUTION	Always centrifuge tubes before opening. Do not mix by vortex or pipetting.





It is not recommended to reconstitute to a concentration less than 100µg/ml.

Dissolve the lyophilized protein in distilled water.

Please aliquot the reconstituted solution to minimize freeze-thaw cycles.

SHIPPING

The product is shipped at ambient temperature.

Upon receipt, store it immediately at the temperature listed below.

STORAGE

Lyophilized protein should be stored at $< -20^{\circ}\text{C}$, though stable at room temperature for 3 weeks.

Reconstituted protein solution can be stored at $4-7^{\circ}\text{C}$ for 2-7 days.

Aliquots of reconstituted samples are stable at $< -20^{\circ}\text{C}$ for 3 months.

BACKGROUND

Proenkephalin-B(PDYN), belongs to the opioid neuropeptide precursor family. The N-terminal domain contains 6 conserved cysteines thought to be involved in disulfide bonding and/or processing. Leu-enkephalins, which is a type of Proenkephalin-B, compete with and mimic the effects of opiate drugs. They play a role in a number of physiologic functions, including pain perception and responses to stress.






ELK Biotechnology

Dynorphin peptides differentially regulate the kappa opioid receptor. Dynorphin A has a typical opioid activity, it is 700 times more potent than Leu-enkephalin.

SDS-PAGE

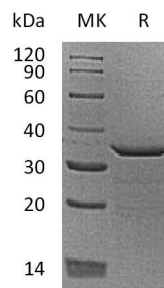

+86-27-59760950

ELKbio@ELKbiotech.com

www.elkbiotech.com

23-2, No.388 Gaoxin 2nd Road, Wuhan East Lake Hi-tech Development Zone, Hubei, P.R.C